



COLLEGE PARK

GUARDIAN™

friction-brake knee



EASY ADJUSTMENTS | USER-FRIENDLY LOCK | STANCE FLEXION




This lightweight, mechanical knee features a friction brake for stance control, stance flexion, and extension assist. It's best suited as a rehabilitation tool for new patients or permanent solution for low-impact users. With a newly designed internal locking mechanism and adjustable toe load sensitivity, patients may find it much easier to lock and unlock. To enhance the sleek look and efficient function, our engineers removed the levers from either side of the knee. All adjustments are now streamlined and internalized!

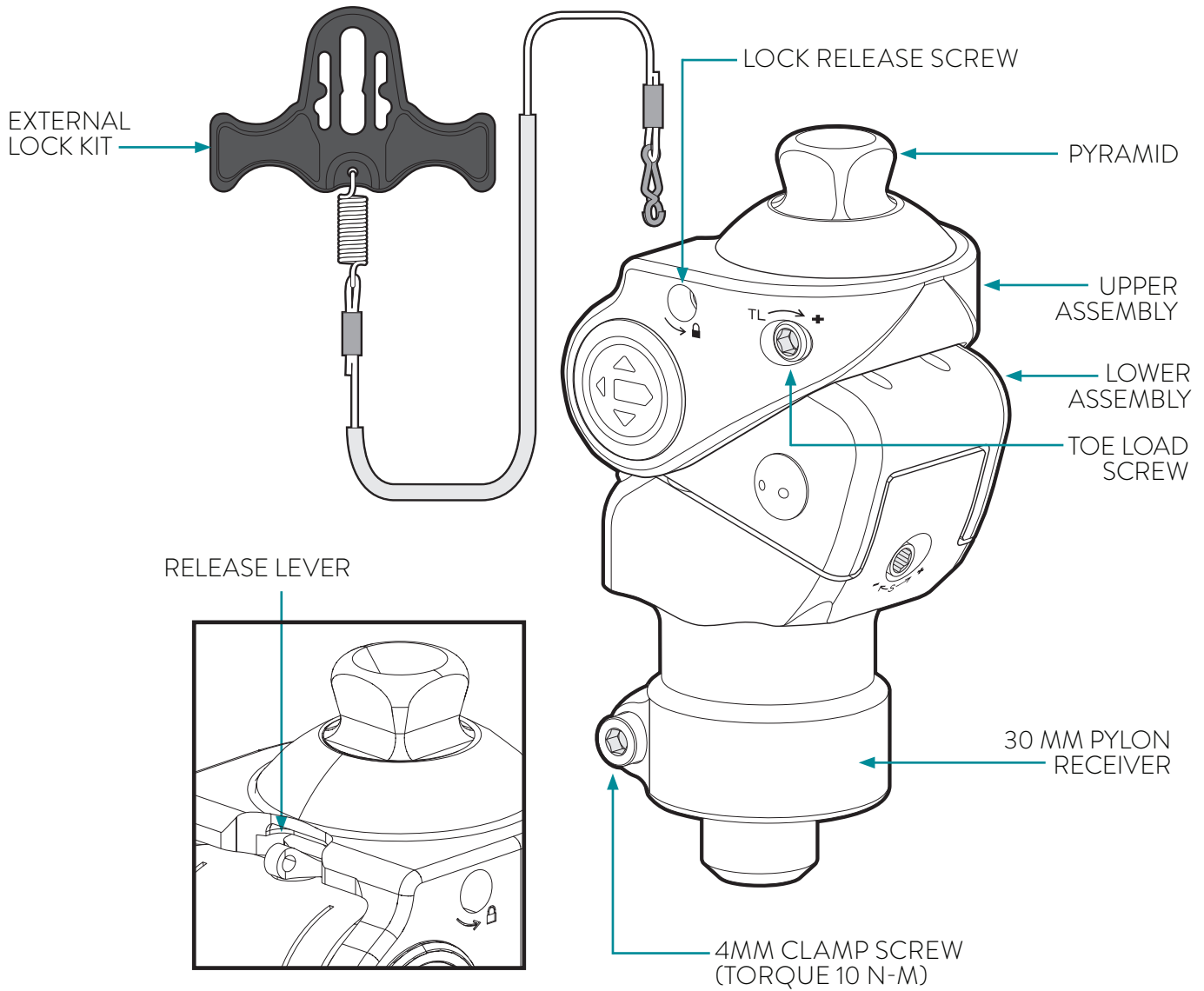
GUARDIAN PART NUMBER

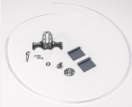

GFKA 120

TECHNICAL SPECIFICATIONS

MATERIAL	Aluminum
DISTAL ATTACHMENT	30 mm Pylon Receiver
FLEXION RANGE	145°
WEIGHT	533 g
WEIGHT LIMIT	275 lbs
WARRANTY	2 years
DOMES TO KNEE CENTER	1.02 in
OVERALL HEIGHT	5.2 in
L-CODES*	L5812, L5845, L5850, L5925
WATER CAPABILITY	
 Weatherproof	

*The listing of codes with these products should not be construed as a guarantee for coverage or payment. Ultimate responsibility for the coding of services/products rests with the individual practitioner.



HARDWARE/TOOLS	PART NUMBER	DESCRIPTION
	CP HK 4	Hex Key, 4mm
	CP KELK	Knee External Lock Kit (included with knee)
	CP LK	Lanyard Kit (included with knee)