

INTEGRATED SPRING TECHNOLOGY | WEDGES FOR FINE TUNING | LOWER PROFILE



The Velocity foot's coupled toe springs work together to help provide a progressive, dynamic response. This foot was designed to enhance smooth rollover for comfortable walking and high impact activities. The Velocity foot combines high functionality for thrill seekers and minimal maintenance to keep them on the go.

INCLUDED

Endo foot includes: foot shell, cosmetic attachment plate (CAP), CPI sock and wedge kit

VELOCITY PART NUMBER

VL	SIDE	SHELL COLOR*	M EN	SIZE	WIDTH†	D	FIRMNESS CATEGORY
	L OR R	C, T, B OR J		25-30	N OR W		1-6

*(C)aucasian, (T)an, (B)rown or (J)et Black. †(N)arrow or (W)ide. Example part number: VL LCMEN 27W D6




ENVIROHELL PART NUMBER

V S	SIDE	SHELL COLOR**	SIZE	WIDTH
	L OR R	C, T, B OR J	25-30	N OR W


** (C)aucasian, (T)an, (B)rown or (J)et Black.



TECHNICAL SPECIFICATIONS

FUNCTION	iST®
SIZES	25-30 cm
MOUNTING	Endo (IP)
WEIGHT LIMIT	250 lbs (25-26 cm) 275 lbs (27-30 cm)
WARRANTY	3 years
BUILD HEIGHT	4.7 in (25-26 cm) 5.2 in (27-30 cm)
WEIGHT*	631 g
L-CODE** L5981, L5986	IMPACT LEVEL 
WATER CAPABILITY  Weatherproof	ENVIROHELL***  narrow & wide

*26 cm with Enviroshell. **The listing of codes with these products should not be construed as a guarantee for coverage or payment. Ultimate responsibility for the coding of services/products rests with the individual practitioner. ***Caucasian, Tan, Brown, or Jet Black.

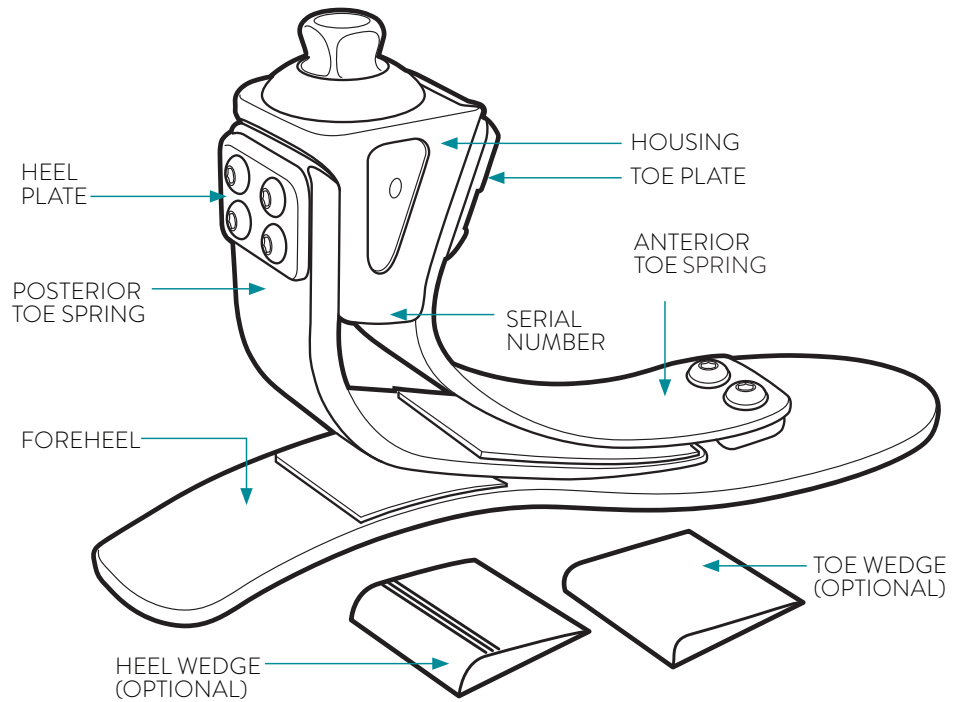

VELOCITY FIRMNESS CATEGORIES

PATIENT WEIGHT	0-100lbs	101-120	121-140	141-160	161-180	181-200	201-220	221-250	251-275
SIZE	25-26cm								
LOW - MODERATE IMPACT	1	1	1	2	2	3	4	5	
HIGH IMPACT	1	2	2	3	3	4	5	6	
SIZE	27-30cm								
LOW - MODERATE IMPACT	1	1	1	1	2	2	3	4	5
HIGH IMPACT	1	2	2	2	3	3	4	5	6

HARDWARE	PART NUMBER	DESCRIPTION
	HZ HWK 2125	Horizon Wedge Kit, 21-25cm
	HZ HWK 2630	Horizon Wedge Kit, 26-30cm
	CP SSET	Super Set Screws (8pc) M8 x 14mm

TOOLS	PART NUMBER	DESCRIPTION
	TFH 800	Foot Horn, Stainless Steel
	CP HK 4	Hex Key, 4mm
	CA GEL	Cyanoacrylate Adhesive Gel (CA glue)

COSMETIC ATTACHMENT PLATE (CAP)	PART NUMBER (LEFT/RIGHT)	DESCRIPTION	ADDITIONAL INFORMATION
	V CLΔ 2526 V CRΔ 2526	Cosmetic Attachment Plate (left or right), 25-26cm	Plastic ring snaps into foot shell for easy cosmetic finishing.
	V CLΔ 2730 V CRΔ 2730	Cosmetic Attachment Plate (left or right), 27-30cm	

Δ Indicate Color. (C)aucasian, (B)rown or (T)an.

SHELLTREAD	PART NUMBER (LEFT/RIGHT)	DESCRIPTION	ADDITIONAL INFORMATION
	V TLN xx V TRN xx	Shelltread (left or right) Neutral (indicate size, width)	Increases durability and prevents premature wear of foot shell.
	V TLJ xx V TRJ xx	Shelltread (left or right) Jet Black (indicate size, width)	

SOCKS	PART NUMBER	DESCRIPTION	ADDITIONAL INFORMATION
	CPI SS SM	CPI Sock, Small, 21-25cm	Highly durable sock that acts as a noise barrier and protects the foot from debris.
	CPI SS LG	CPI Sock, Large, 26-30cm	
	CPI SSK SM	CPI Sock, Small, 21-25cm, 6-pack	
	CPI SSK LG	CPI Sock, Large, 26-30cm, 6-pack	