

TECHNICAL SPECIFICATIONS



FUNCTION	MOUNTING	WEIGHT LIMIT		CLEARANCE		IMPACT LEVEL	WATER CAPABILITY
Tri-Axial	Endo (IP) Exo (ALX)	21-24 cm	250 lbs (113 kg)	Endo (IP)	2.5-2.7 in (6.4-6.9 cm)	L M H K3 K4	WEATHERPROOF
SIZES	WEIGHT*	WARRANTY	L CODES**				
21-30 cm	585 g	3 years	L5981 + L5986		ENVIROSHELL®*** CTBD narrow & wide		

VENTURE PART NUMBER

SIDE	SHELL COLOR**	HOUSING†	SIZE	WIDTH‡	FIRMNESS CATEGORY		
V	L or R	C, T, B or J	EN or ALX	21-30	N or W	D	1-5

*** (C)caucasian, (T)an, (B)rown or (J)et Black. † (EN)do, Exo (ALX). ‡ (N)arrow or (W)ide.
Example part number: V LCEN 24W D4



FUNCTION	MOUNTING	WEIGHT LIMIT		CLEARANCE		IMPACT LEVEL	WATER CAPABILITY
Tri-Axial	Endo (IP) Exo (ALX)	Endo (IP)	330 lbs (150 kg)	Endo (IP)	3.5 in (8.7 cm)	L M H K3 K4	WEATHERPROOF
SIZES	WEIGHT*	WARRANTY	L CODES**				
22-31 cm	678 g	3 years	L5979 PDAC APPROVED		ENVIROSHELL®*** CTBD narrow & wide		

TRUSTSTEP PART NUMBER

SIDE	SHELL COLOR***	HOUSING	SIZE	WIDTH‡	FIRMNESS CATEGORY		
TS	L or R	C, T, B or J	EN or ALX	22-31	N or W	D	1-9

*** (C)caucasian, (T)an, (B)rown or (J)et Black. ** (EN)do, Exo (ALX). † (N)arrow or (W)ide. Wide only available in 24-31cm. Example part number: TS LCEN 24N D6



FUNCTION	MOUNTING	WEIGHT LIMIT		CLEARANCE		IMPACT LEVEL	WATER CAPABILITY
Tri-Axial	Endo (IP)	220 lbs (100 kg)		2.4 in (6.1 cm)		L M K3	WEATHERPROOF
SIZES	WEIGHT*	WARRANTY	L CODES**				
21-30 cm	529 g	2 years	L5979 PDAC APPROVED		ENVIROSHELL®*** CTBD standard narrow & wide		

TRIBUTE PART NUMBER

SIDE	SHELL COLOR***	SANDAL TOE†	SIZE	WIDTH‡	FIRMNESS CATEGORY			
TB	L or R	C, T, B or J	EN	S	21-30	N or W	D	1-4

*** (C)caucasian, (T)an, (B)rown or (J)et Black. † (S)sandal. Sandal toe not available in Jet Black. Leave blank for regular toe. ‡ (N)arrow or (W)ide. Wide not available for Sandal Toe. Example part number: TB LCEN 24N D3

*26 cm model, including Envirosheel. **The listing of codes with these products should not be construed as a guarantee for coverage or payment. Ultimate responsibility for the coding of services/products rests with the individual practitioner. *** Caucasian, Tan, Brown, or Jet Black. † Low impact only.

ALSO AVAILABLE!


SERVICE FREE



TRI-AXIAL TECHNOLOGY


Tri-axial feet offer movement across all three planes: the sagittal, coronal, and transverse. Our journey began in 1988 when we developed the first prototype of the anatomically correct prosthetic foot, the Truststep. Drawing on decades of expertise in tri-axial technology, our team designed more solutions that empower users to navigate uneven terrain with ease.

THREE PLANES OF MOVEMENT

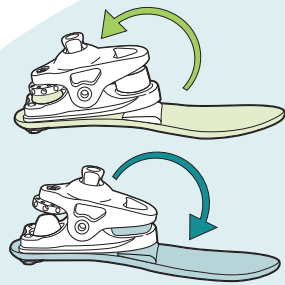


PLANTARFLEXION/ DORSIFLEXION

These feet were created to maximize ground contact. In addition to movement from the bumpers in the ankle unit, tri-axial feet achieve plantar and dorsiflexion through springs.

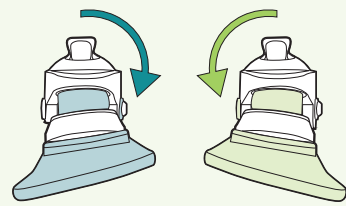


INVERSION/ EVERSION

From hiking a trail to exploring the city, tri-axial feet were created to accommodate bumps in the road with excellent inversion and eversion.

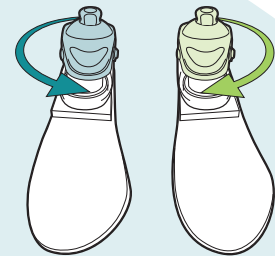

TRANSVERSE ROTATION

Users may benefit from increased comfort as they twist and turn throughout the day.

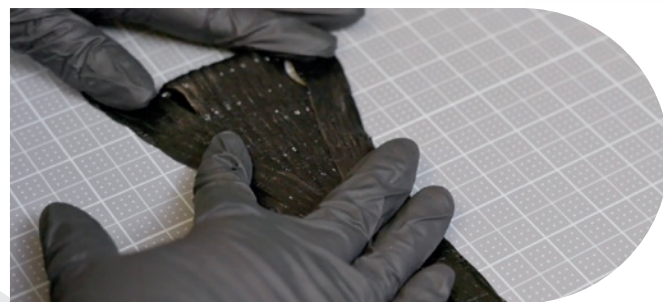

MATERIALS MATTER: INTELLIWEAVE

College Park's Intelliweave® fiberglass composite technology aims to enhance flexibility and comfort. These fibers are woven by hand in a 3D pattern for strength and precision gait matching. Intelliweave is responsive to activity, which can make each step feel more natural while reducing impact on the user's limb and socket.



BUMPER FLEXIBILITY

College Park offers a variety of urethane bumpers and bushings for optimal customization. These long-lasting bumpers are made of aircraft grade polymers for enhanced durability.



L5981, L5986**

THE ULTIMATE FUSION

With a full-length toe lever, the Venture foot was designed to provide a higher frequency dynamic response for more active users. Highly engineered bumpers and flexible Intelliweave® composites work together for a progressive response. The tri-axial, custom-built foot includes our exclusive Stride Control™ feature for easy fine-tuning without disassembly.



CUSTOMIZATION



"The Venture foot is very durable and long-lasting. I've had this one for quite a few years now. I can walk, run, skateboard, snowboard—it's really all encompassing. I wouldn't be able to do what I do without it."

— Oscar, Venture Foot User



L5979**
PDAC APPROVED

IDEAL CHOICE FOR EVERYDAY LIFE

Original and unrivaled, the Truststep foot sets the standard for comfort and natural-feeling mobility. It was one of the first to provide correct movement in all planes—allowing for nearly the same vertical motion, rotation, and stability as an anatomical foot. With a full range of options available, the versatile design allows complete customization.



6.8 MILLIMETERS OF
VERTICAL COMPLIANCE



POLYCENTRIC



"After I put my legs on in the morning, I don't take them off until I go to bed in the evening. I've been told that I have a great gait. As a longtime amputee, I've worn many legs, and it hasn't always been like that. I started wearing my Truststeps in 2001 and haven't looked back."

— Reggie, Truststep Foot User



L5979**
PDAC APPROVED

SOCIAL MOBILITY

The Tribute foot delivers dependable performance in a value package. Lightweight and low-profile, it uses a full-length toe lever and tri-axial motion to help achieve firm footing on uneven ground. This foot is an ideal solution for low to moderate impact users.



CUSTOMIZATION



**The listing of codes with these products should not be construed as a guarantee for coverage or payment.