


STERLING

Gel Liner

The College Park Sterling is a seamless TPE gel liner designed to conform to the patient's limb to provide optimal comfort. The Thermoplastic Elastomer (TPE) gel can accommodate limb irregularities. This durable liner can be easy to don and doff with an optimal TPE surface. It features a soft, extra thick distal pad which may improve user comfort, protect sensitive distal ends, and conform to irregular distal shapes.

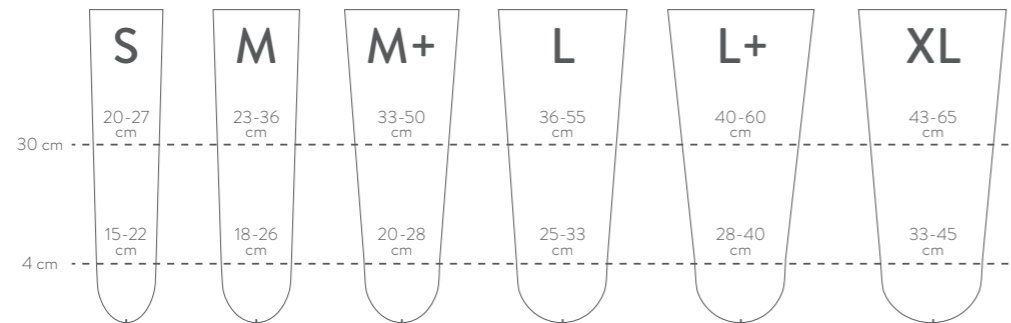

STERLING LINER PART NUMBER

SLK **G** **3,6,OR9** **LORC** **U** **S-XL** **1OR2**

*See charts above. Pin (L)ock or (C)ushion.

TECHNICAL SPECIFICATIONS

MATERIAL	TYPES	PROFILES
Thermoplastic Elastomer (TPE) Gel	Cushion or Locking	Uniform
SIZES	THICKNESSES	OVERALL LENGTH
S, M, M+, L, L+, XL	3 mm, 6 mm, 9 mm	16.5 in (42.0 cm)



STERLING

Gel Sleeve

The new Sterling sleeve features Thermoplastic Elastomer (TPE) Gel, designed to use body heat and mold to a user's limb for a near-custom fit. It features a seamless textile to help prevent snagging, with a thicker middle section for durability.



WHY TPE?

The Thermoplastic Elastomer (TPE) gel uses body heat to mold to a user's limb for a custom fit. Designed to accommodate limb irregularities and enhance comfort, TPE could make the Sterling an ideal solution for low to moderate impact users.

TO DETERMINE THE CORRECT SIZE:

Measure the circumference of the residual limb 15 cm above mid-patella and choose the corresponding device size. See the sizing chart below.

SIZE	MEASUREMENT (CM)	PART NUMBER	DESCRIPTION
35	35-41 cm	SGS 3541	Sterling Suspension Sleeve, Size 35-41 cm
41	41-49 cm	SGS 4149	Sterling Suspension Sleeve, Size 41-49 cm
49	49-55 cm	SGS 4955	Sterling Suspension Sleeve, Size 49-55 cm
55	55-62 cm	SGS 5562	Sterling Suspension Sleeve, Size 55-62 cm
62	62-76 cm	SGS 6276	Sterling Suspension Sleeve, Size 62-76 cm

