

# TEMPO® iVT

integrated vacuum technology  
intelliweave composite  
service-free design



## TECH SPECS

### FUNCTION

K2 Vacuum

### BUILD HEIGHT

2.5 in

### MOUNTING

Endo (IP)

### WEIGHT LIMIT

21-24 cm | 300 lbs

25-30 cm | 330 lbs

### L-CODE†

L5972

L5781

### WATER CAPABILITY



WEATHERPROOF

### SIZES

21-30 cm

### WEIGHT\*

626 g

### IMPACT LVL



### WARRANTY

2 years

### ENVIROSHIELD†



standard



For people seeking vacuum benefits with a higher weight limit, College Park offers the Tempo iVT! This low-profile, lightweight foot features College Park's Integrated Vacuum Technology (iVT), which is designed to pull up to 20 inHg in fewer steps.\* This allows patients to efficiently achieve optimal vacuum. The service-free design comes with a full length toe lever, providing stored energy for a natural response and extra ground contact for stability. The Tempo iVT is designed to enhance user confidence, comfort, and more.

\*Results may vary based on impact and activity. Maximum attainable vacuum during full dorsiflexion ranges from 18-22 inHg.

## TEMPO iVT PART NUMBER

	SIDE	SHELL COLOR*		SIZE		FIRMNESS CATEGORY
<b>TO IVT</b>	<b>L OR R</b>	<b>C, T OR B</b>	<b>EN</b>	<b>21-30</b>	<b>D</b>	<b>1-5</b>

\*(C)aucasian, (T)an or (B)rown. Shell width is standard. Example part number: TO IVT LCEN 24 D3

## TEMPO iVT FIRMNESS CATEGORIES



PATIENT WEIGHT	0-140lbs	141-180	181-220	221-300	301-330
SIZE	21-30cm				25-30cm
LOW IMPACT	1	2	3	4	5


## INCLUDED

Endo foot with integrated vacuum pump, includes: bonded sandal toe foot shell, CPI sock, and choice of AK or BK tube kit.


## TUBE KITS

(Included with foot order, selection required)


PART NO	DESCRIPTION	
 IVT AK TBK	iVT AK Tube Kit	Includes: AK tube, tube retainer, small elastic band, large elastic band, syringe, hose T-connector.
 IVT BK TBK	iVT BK Tube Kit	Includes: BK tube, tube retainer, syringe, hose T-connector.

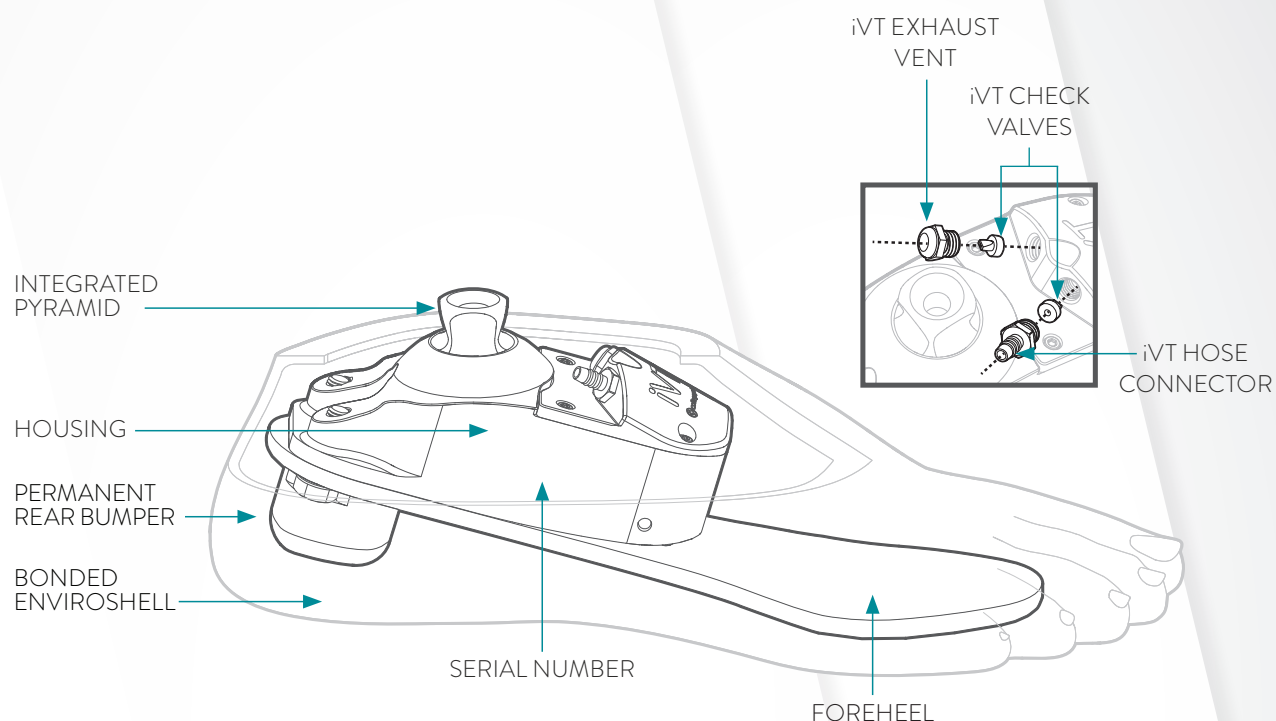

## ACCESSORIES

	PART NUMBER	DESCRIPTION	ADDITIONAL INFORMATION
	IVT VAK	iVT Valve Assembly Kit	Tri-function socket valve. Includes: iVT valve, 1 oz RTV silicone, iVT valve air filter pack (10pc), cap key, socket insert, adhesive pad, screw.

## SERVICE KITS (for servicing replaceable portions of vacuum pump)

	PART NUMBER	DESCRIPTION	ADDITIONAL INFORMATION
	CTV SVK 2124	CT Vacuum Service Kit, 21-24cm	All parts needed to perform full service on replaceable portions of vacuum pump. Includes: vacuum bladder, (2) return spring assemblies, (2) check valves, hose connector, exhaust vent, M3 x 0.50mm air bladder fastener, (4) M2.5 x 6.5 housing cap fasteners, BK tube, T9 Torx wrench, 5/16 wrench.
	CTV SVK 2526	CT Vacuum Service Kit, 25-26cm	
	CTV SVK 2730	CT Vacuum Service Kit, 27-30cm	





TEMPO iVT


## REPLACEMENT PARTS

	PART NUMBER	DESCRIPTION	ADDITIONAL INFORMATION
	IVT AK VK	iVT AK Valve Kit	AK Tube Kit + iVT Valve Assembly Kit. Includes: AK tube, tube retainer, small elastic band, large elastic band, syringe, hose T-connector, iVT valve, 1 oz RTV silicone, iVT valve air filter pack (10pc), cap key, socket insert, adhesive pad, screw.
	IVT BK VK	iVT BK Valve Kit	BK Tube Kit + iVT Valve Assembly Kit. Includes: BK tube, tube retainer, syringe, hose T-connector, iVT valve, 1 oz RTV silicone, iVT valve air filter pack (10pc), cap key, socket insert, adhesive pad, screw.
	IVT AK TB	iVT AK Tube	
	IVT BK TB	iVT BK Tube	
	CTV VB 2124	CT Vacuum Bladder, 21-24cm	
	CTV VB 2526	CT Vacuum Bladder, 25-26cm	
	CTV VB 2730	CT Vacuum Bladder, 27-30cm	
	CTV RSK 2124	CTV Return Spring Assembly Kit (2 pc) 21-24cm	
	CTV RSK 2526	CTV Return Spring Assembly Kit (2 pc) 25-26cm	
	CTV RSK 2730	CTV Return Spring Assembly Kit (2 pc) 27-30cm	
	IVT CKV	iVT Check Valve	
	IVT HC	iVT Hose Connector	
	IVT EXV	iVT Exhaust Vent	
	CTV BFK	CTV Bladder Fastener Kit	
	CTV HCFK	CT Vacuum Housing Cap Fastener Kit (4pc)	
	IVT VIC	iVT Valve Inner Core	iVT Valve inner core only.
	IVT VFP	iVT Valve Air Filter Pack (10pc)	Replaceable particulate filters for iVT Valve. Includes: (10) filters.

## HARDWARE AND TOOLS

	PART NUMBER	DESCRIPTION
	CP SSET	Super Set Screws (8pc) M8 x 14mm
	CP HK 4	Stride Control Key, 4mm

## SHELLTREAD

	PART NUMBER (LEFT/RIGHT)	DESCRIPTION	ADDITIONAL INFORMATION
	CS TLN x      CS TRN x	Shelltread (left or right), Neutral (indicate size)	Increases durability and prevents premature wear of foot shell.