

TECHNICAL SPECIFICATIONS



FUNCTION	MOUNTING	WEIGHT LIMIT	CLEARANCE		IMPACT LEVEL	WATER CAPABILITY
Adjustable Heel	Endo (IP)	220 lbs (21-24 cm)	IP	3.3-3.9 in (8.5-9.9 cm)		 WEATHERPROOF
	D.Pylon (DP)	250 lbs (25-28 cm)	4-Hole	7.8-8.0 in (19.8-20.2 cm)		
SIZES	WEIGHT*	WARRANTY	L CODES**		ENVIROSHELL***	
21-28 cm	(IP)	608 g	L5981, L5986, L5990			
	(DP) w/pylon	761 g				
		2 years				

ACCENT PART NUMBER

AC	SIDE	SHELL COLOR***	MOUNTING	BUTTON PLACEMENT	SIZE	WIDTH‡	FIRMNESS CATEGORY
	L or R	C or B	IP, DF, DM or D4	M or L	21-28	N or W	D 1-3

*(C)aucasian or (B)rown. †(IP) Integrated Pyramid, (DF) Dynamic Pylon Female, (DM) Dynamic Pylon Male or (D4) Dynamic Pylon 4-hole. **(M)edial placement or (L)ateral placement. ‡(N)arrow or (W)ide. Wide only available in 23-26cm. Example part number: AC LCIP 24N D2



FUNCTION	MOUNTING	WEIGHT LIMIT	CLEARANCE	IMPACT LEVEL	WATER CAPABILITY
Stubbies	Endo (IP)	275 lbs (125 kg)	2.3 in (5.8 cm)		 WATERPROOF
SIZES	WEIGHT	WARRANTY	L CODES**		
One size	422 g	2 years	L5972, L5986		

SIDEKICKS PART NUMBER

SK	AK	QUANTITY*	FIRMNESS CATEGORY
		1 or 2	D 1-3

*1 for single, 2 for pair. Example part number: SK AK2 D2

*25 cm model, including Enviroshell. **The listing of codes with these products should not be construed as a guarantee for coverage or payment. Ultimate responsibility for the coding of services/products rests with the individual practitioner. ***Caucasian or Brown.

ALSO AVAILABLE!



ENDO COMPONENTS



SPECIALTY: ADJUSTABLE HEEL HEIGHT



L5981 + L5986 + L5990..

DESIGNED FOR FASHION. ENGINEERED FOR LIFE.

The Accent foot provides 2” (5 cm) of heel height adjustment through a simple push of a button. It comes with a contoured, soft shin fairing that locks into the shell and prevents cosmetic buckling when the foot is adjusted at various heel heights. It is a service-free option for individuals who desire a cosmetically appealing foot.



SERVICE-FREE
DESIGN




COSMETIC
OPTIONS

MORE RANGE THAN EVER

The Accent has 13 increments of adjustment that lock securely into place. With an easy-to-adjust button located in the ankle, the foot accommodates a wide range of shoes including flats, cowboy boots, and high-heels.



ADJUSTABLE HEEL HEIGHT THROUGH THE SIMPLE PUSH OF A BUTTON!

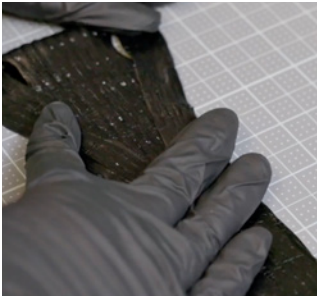

MATERIALS MATTER: INTELLIWEAVE

College Park’s Intelliweave® composite technology aims to enhance flexibility and comfort. These fibers are woven by hand in a 3D pattern for strength and precision gait matching. Intelliweave is responsive to activity, which can make each step feel more natural while reducing impact on the user’s limb and socket. The Accent foot uses Intelliweave composites to help maximize performance in everyday life.



“I’m pulling out all of my boots from the back of my closet! After almost a year of only wearing sneakers, flat shoes, zero heel, I can now adjust the angle of my foot up to 2 inches. It’s pretty exciting.”

- Karen, Accent User

SPECIALTY: FLEXIBLE STUBBIES



L5972 + L5986..

UP, UP, AND AWAY

The Sidekicks® are adjustable stubbie feet for foreshortened prostheses. College Park engineers applied 15 years of knowledge gained from the Venture® ankle technology to develop the Sidekicks. With anatomical ankle motion and a flexible platform, they were designed to be a vast improvement over other options for bilateral, transfemoral amputees. Gait-matched and fully waterproof, this foot can flex and plant in real-world environments for more balance and stability.



STRIDE
CONTROL



VERSATILE
FUNCTION



GAIT
MATCHED

PLANTARFLEXION



DORSIFLEXION



INVERSION/EVERSION



TRANSVERSE ROTATION

ANKLE ROTATION

The Sidekicks feet were designed to mimic anatomical ankle movement, which may help reduce socket forces and eliminate unwanted torque.



REHABILITATION TO RECREATION

With natural ankle movement on a flexible platform, Sidekicks were created to provide comfort through rehabilitation or daily life. They’re fully waterproof and work well for many recreational activities!



“While competing in New Zealand, I ran them through four airports, three cities, and two ski resorts, through sand, snow, water, mud, and more. They held up to all conditions.”

- Josh, Sidekicks User

**The listing of codes with these products should not be construed as a guarantee for coverage or payment.