

# VOYAGER®

vertical shock absorption  
smooth rotation  
material fusion



The Voyager was designed to achieve ultimate comfort and rotation for any journey. With responsive springs, vertical shock absorption, and built-in torsion, this high-performance foot can feel more natural to the user. Its superior range of motion accommodates varying activity levels. The sleek Voyager uses a hybrid material approach—blending flexible Intelliweave® composites with durable carbon fiber to maximize function.

K3/4



## TECH SPECS

### FUNCTION

Vertical Shock +  
Torsion

### BUILD HEIGHT

7.5 in

### MOUNTING

Endo (IP)

### WEIGHT LIMIT

21-25 cm | 275 lbs

26-30 cm | 330 lbs

### L-CODE†

L5987, L5984

### SIZES

21-30 cm

### WEIGHT\*

1,302 g

### IMPACT LVL



### WARRANTY

3 years

### ENVIROSHELL†



standard

\*26 cm with Enviroshell. ‡ The listing of codes with these products should not be construed as a guarantee for coverage or payment. Ultimate responsibility for the coding of services/products rests with the individual practitioner. †Caucasian, Tan, Brown and Jet Black.

## VOYAGER PART NUMBER

	SIDE	SHELL COLOR*		SIZE		FIRMNESS CATEGORY
<b>VY</b>	<b>L OR R</b>	<b>C, T, B OR J</b>	<b>EN</b>	<b>21-30</b>	<b>D</b>	<b>1-5</b>

\*(C)aucasian, (T)an, (B)rown or (J)et Black. ‡(N)arrow or (W)ide. Example part number: VY LBEN 21 D1

## VOYAGER FIRMNESS CATEGORIES

PATIENT WEIGHT	0-140lbs	141-180	181-220	221-275	276-330
SIZE	21-25cm				
MODERATE-HIGH IMPACT	1	2	3	4	
SIZE	26-30cm				
MODERATE-HIGH IMPACT	1	2	3	4	5

## ENVIROSHELL PART NUMBER

	SIDE	SHELL COLOR**	SIZE
<b>HZ S</b>	<b>L OR R</b>	<b>C, T, B, OR J</b>	<b>21-30</b>

\*\* (C)aucasian, (T)an, (B)rown, or (J)et Black




Enviroshell included with every foot order



## INCLUDED

Endo foot includes: foot shell, CPI sock and wedge kit

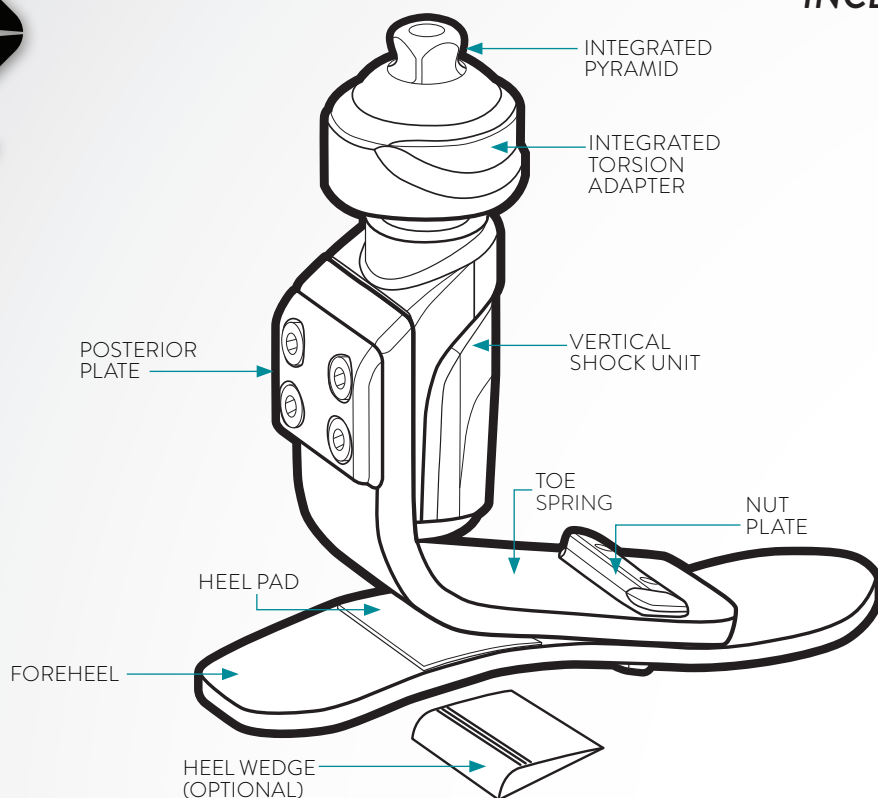

## HARDWARE

	PART NUMBER	DESCRIPTION
	VY HWK 2125	Voyager Wedge Kit, 21-25
	VY HWK 2630	Voyager Wedge Kit, 26-30cm
	CP SSET	Super Set Screws (8pc) M8 x 14mm

## TOOLS

	PART NUMBER	DESCRIPTION
	TFH 800	Foot Horn, Stainless Steel
	CP HK 4	Hex Key, 4mm
	CA 409	Cyanoacrylate Adhesive, 409 Gel (CA glue)

## SHELLTREAD

	PART NUMBER (LEFT/RIGHT)	DESCRIPTION	ADDITIONAL INFORMATION
	HZ TLN xx      HZ TRN xx	Shelltread (left or right) Neutral (indicate size, width)	Increases durability and prevents premature wear of foot shell.
	HZ TLJ xx      HZ TRJ xx	Shelltread (left or right) Jet Black (indicate size, width)	

## SOCKS

	PART NUMBER	DESCRIPTION	ADDITIONAL INFORMATION
	HZ SS 2125	CPI Sock, 21-25cm	Highly durable sock that acts as a noise barrier and protects the foot from debris.
	HZ SS 2630	CPI Sock, 26-30cm	
	HZ SSK 2125	CPI Sock, 6-pack, 21-25cm	
	HZ SSK 2630	CPI Sock, 6-pack, 26-30cm	