


# SIDEKICKS®

tri-axial function  
adjustable gait match  
molded tread

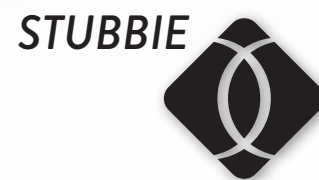

## TECH SPECS

### FUNCTION

Stubbies

### BUILD HGT

2.3 in

### MOUNTING

Endo (IP)

### WEIGHT LIMIT

275 lbs

### L-CODE†

L5972  
L5986

### SIZES

One size

### WEIGHT

422 g

### IMPACT LVL



### WARRANTY

2 years

### WATER CAPABILITY



†The listing of codes with these products should not be construed as a guarantee for coverage or payment. Ultimate responsibility for the coding of services/products rests with the individual practitioner.



Sidekicks are the first adjustable stubbie feet for foreshortened prostheses. College Park engineers applied 15 years of knowledge gained from the Venture ankle technology to develop the Sidekicks. With natural anatomical ankle motion and a flexible platform for ground compliance, they are a vast improvement over current options for bilateral, transfemoral amputees.

## SIDEKICKS PART NUMBER

**SK AK** **QUANTITY\*** **1 OR 2** **D** **FIRMNESS CATEGORY** **1-3**

\*1 for single, 2 for pair. Example part number: SK AK2 D2

## SIDEKICKS FIRMNESS CATEGORIES


PATIENT WEIGHT	0-120lbs	121-220	221-275
SIZE	ONE SIZE		
LOW IMPACT	1	2	3

## INCLUDED

Endo foot includes:  
Sidekicks bumper kit  
and pin guide tool

## GAIT MATCHING FIRMNESS CATEGORY



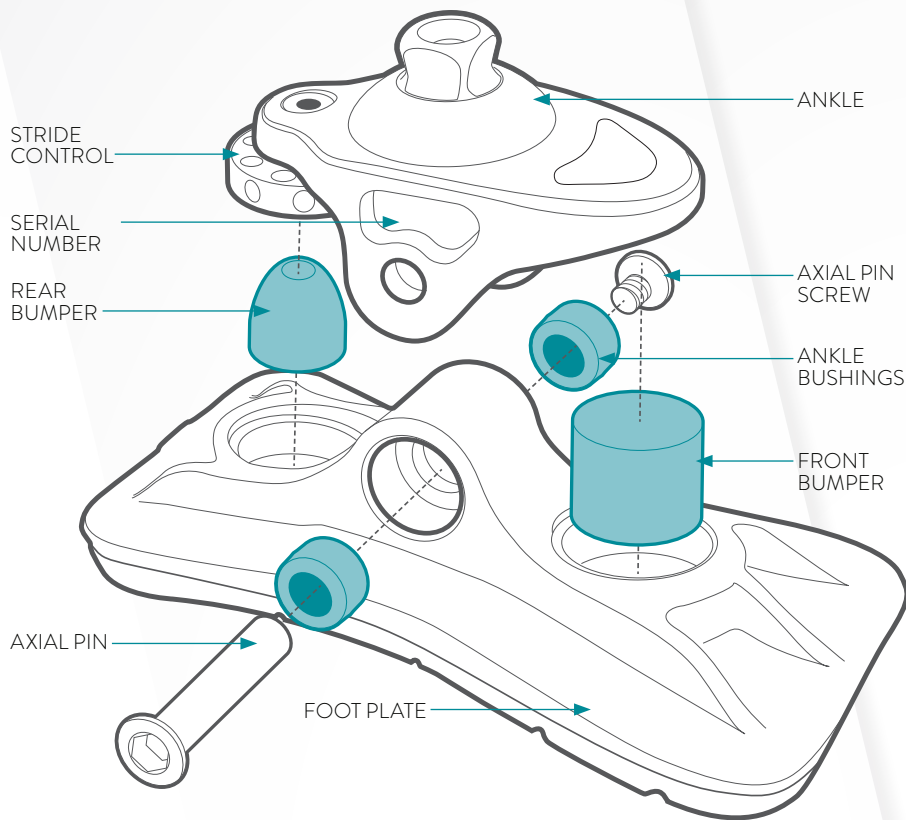


COMPONENT	1	2	3
Front Bumper	S	M	F
Ankle Bushings	M	M	M
Rear Bumper	F	F	F



PART NO	DESCRIPTION
SK FP	Sidekicks Foot Plate



## SOFT COMPONENTS

PART NO	DESCRIPTION
 V FBS A	Front Bumper, Soft, Size A
 V FBM A	Front Bumper, Medium, Size A
 V FBF A	Front Bumper, Firm, Size A
 V RBS B	Rear Bumper, Soft, Size B
 V RBM B	Rear Bumper, Medium, Size B

PART NO	DESCRIPTION
 V RBF B	Rear Bumper, Firm, Size B
 V RBX B	Rear Bumper, Extra Firm, Size B
 V ABS 2530	Ankle Bushing, Soft, 25-30cm
 V ABM 2530	Ankle Bushing, Medium, 25-30cm
 V ABF 2530	Ankle Bushing, Firm, 25-30cm

PART NO	DESCRIPTION
 SK BK 1	Sidekicks Bumper Kit, Single (See page 52 for contents)
 SK BK 2	Sidekicks Bumper Kit, Pair (See page 52 for contents)


FIRMNESS  
COLOR  
CODE:

 SOFT  
 MEDIUM  
 FIRM

## HARDWARE

PART NUMBER	DESCRIPTION
 CP AP L	Axial Pin, Long, 25-30cm
 CP APS	Axial Pin Screw
 CP SSET	Super Set Screws (8pc) M8 x 14mm

## TOOLS

PART NUMBER	DESCRIPTION
 TGP 525	Pin Guide
 CP HK 4	Hex Key, 4mm
 TL LUB	TruLube Lubricant, 0.25 oz.