



PRODUCT CATALOG



Partial Hand
SOLUTIONS, LLC


FULL M-FINGERS




M-Fingers are simple to operate and controlled by the user's own motion, allowing the patient to perform tasks more effectively, more responsively, and with optimized precision. They conform around objects to provide a secure and confident grip for the user. M-Fingers were designed with high strength, injected molded composite materials to deliver the durability and longevity needed for daily use.

FULL M-FINGER SIZING


Ultra-Small M-Finger
Total: 3.19" (80.8 mm)

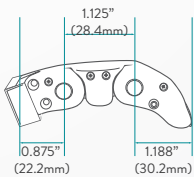


Small M-Finger
Total: 3.63" (91.8 mm)

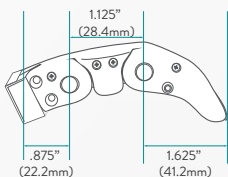


Medium M-Finger
Total: 3.88" (98.3 mm)

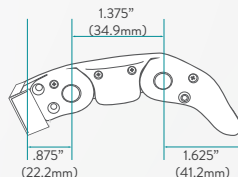

Large M-Finger
Total: 4.25" (107.7 mm)

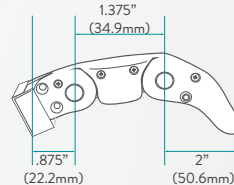

X-Large M-Finger
Total: 4.5" (114 mm)

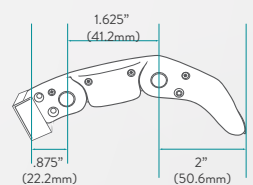

TAN



BROWN



BLACK

M-FINGER PART NUMBER

SIZE

ULTRA-SM, SM, MED, LG, XL

- MF -

COLOR*

TAN, BRN, BLK

* (TAN) Tan, (BRN) Brown, or (BLK) Black. Example part number: ULTRA-SM-MF-BLK

M-THUMB

The M-Thumb maximizes dexterity and functional return for patients with a complete absence of a thumb. It provides a considerable improvement in function due to the multi-positional capabilities and adjustable resistance. The M-Thumb can be rotated and flexed or extended to a position most useful for the user. A fingernail is incorporated into the design, allowing the user to pick up very small objects with ease and precision.



THUMB PART NUMBER

COLOR*

THUMB -

TAN, BRN, BLK

*(TAN) Tan, (BRN) Brown, or (BLK) Black. Example part number: THUMB-BLK

PARTIAL M-FINGERS


The Partial M-Fingers were designed for someone missing a portion of their finger. A custom molded, high strength carbon fiber socket can be constructed to provide optimum comfort and stability. Partial M-Fingers are controlled by motion in the remaining finger so the user has direct control of speed and placement of the fingertip, allowing for precise and rapid grasping of objects.



PARTIAL M-FINGER SIZING


Small Partial M-Finger

Total: 1.65"

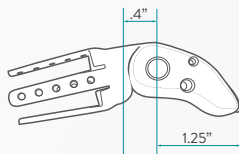


Medium Partial M-Finger

Total: 2.03"

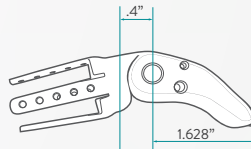

Large Partial M-Finger

Total: 2.4"

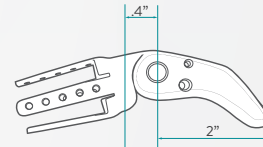

TAN



BROWN



BLACK

PARTIAL M-FINGER PART NUMBER

SIZE

SM, MED, LG

COLOR*

TAN, BRN, BLK

- PART -

* (TAN) Tan, (BRN) Brown, or (BLK) Black. Example part number: SM-PART-BLK


PEDIATRIC M-FINGERS

Developed by Partial Hand Solutions, these fingers are simple to operate and controlled by the child's own movements. This makes the fingers more effective, responsive, and allows for optimized precision. They conform around objects to provide a secure and confident grip.




PEDIATRIC M-FINGER SIZING


Small Partial M-Finger
Total: 2.05" (52.07 mm)

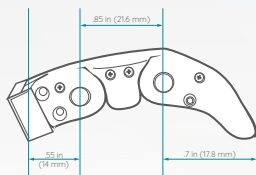

Medium Partial M-Finger
Total: 2.4" (61 mm)

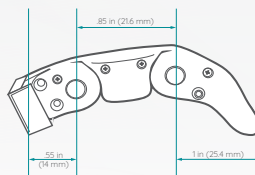

Large Partial M-Finger
Total: 2.65" (67.3 mm)

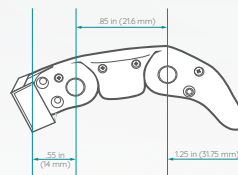



PEDIATRIC M-FINGER SYSTEM PART NUMBER

Pediatric 3 Finger M-Finger System with 3-Finger Mounting Plate, DCM, and BOA mechanism

SIDE

3FS - KIT **L OR R** **- PED**

Example part number: 3FS-KIT-L-PED



PEDIATRIC M-FINGER SYSTEM PART NUMBER

Pediatric 4 Finger M-Finger System with 4-Finger Mounting Plate, DCM, and BOA mechanism

SIDE

4FS - KIT **L OR R** **- PED**

Example part number: 4FS-KIT-L-PED



PEDIATRIC M-FINGER SYSTEM PART NUMBER

5 Finger M-Finger System with 5-Finger Mounting Plate, DCM, and BOA mechanism

SIDE

5FS - KIT **L OR R** **- PED**

Example part number: 5FS-KIT-L-PED

KITS



4 FINGER SYSTEM PART NUMBER

4 Finger M-Finger System with 4-Finger Mounting Plate, DCM, and BOA mechanism

4FS- KIT - COLOR* **TAN, BRN, BLK** **-** SIZE **SM, MED, LG**

* (TAN) Tan, (BRN) Brown, or (BLK) Black. Example part number: SM-PART-BLK



5 FINGER SYSTEM PART NUMBER

4 Finger M-Finger System with 4-Finger Mounting Plate, Thumb, DCM, and BOA mechanism

5FS- KIT - COLOR* **TAN, BRN, BLK** **-** SIZE **SM, MED, LG**

* (TAN) Tan, (BRN) Brown, or (BLK) Black. Example part number: 5FS-KIT-BLK-SM

HOLD-IT

The Hold-It from Partial Hand Solutions allows patients with complete finger amputations or an impaired intact hand the ability to hold tools such as writing implements, eating utensils, and personal care items. The Hold-It easily snaps open and closed, and it can be rotated to accommodate a myriad of different positions. Fabrication is simple and easy. Each Hold-It order includes all the necessary hardware for convenient setup.



HOLD-IT SYSTEM PART NUMBER

HOLD-IT System with Mounting

HOLD-IT


TITAN SERIES




From Partial Hand Solutions, the titanium Titan series fingers are designed for the most active heavy duty user, ideal for those with occupational requirements. They are suitable for someone with a full or partial finger amputation and are available in a variety of sizes. The ratchet design allows the patient to manually position the joints in the most functional position, while providing a secure grasp. Various bases allow for the device to be fit as a thumb with full rotation and adjustable friction.

TITAN FINGER SIZING


TITAN Full 1:

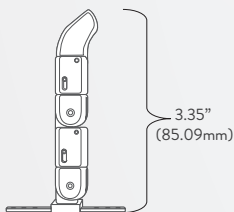


TITAN Full 2:

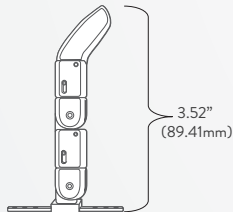


TITAN Full 3:

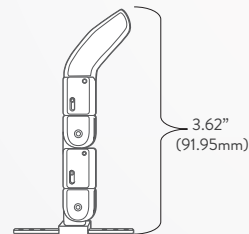


TITAN Partial 1:

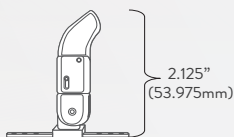


TITAN Partial 2:

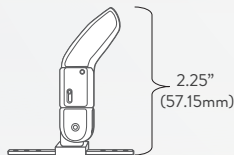


TITAN Partial 3:

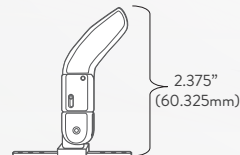


TITAN Thumb 1:

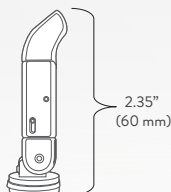

TITAN Thumb 2:

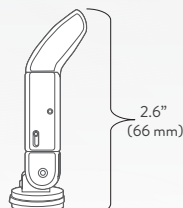

TITAN Thumb 3:

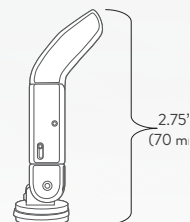



TITAN FINGER PART NUMBER

SIZE[†]

FULL - TITAN

1, 2, 3

[†] (1) Small: length size 3.35", (2) Medium: length size 3.52", or (3) Large: length size 3.62".
Example part number: FULL-TITAN3



TITAN PARTIAL FINGER PART NUMBER

SIZE[†]

PART - TITAN

1, 2, 3

[†] (1) Small: length size 2.25", (2) Medium: length size 2.5", or (3) Large: length size 2.62".
Example part number: PART-TITAN3



TITAN THUMB PART NUMBER

SIZE[†]

THUMB - TITAN

1, 2, 3




[†](1) Small: length size 2.35", (2) Medium: length size 2.6", or (3) Large: length size 2.75".
Example part number: THUMB-TITAN3

ACCESSORIES AND PARTS

MOUNTING PLATES

	PRODUCT NUMBER	DESCRIPTION
	4F-MP	4 Finger Mounting Plate (in-line)
	3F-MP	3 Finger Mounting Plate (in-line)
	2F-MP	2 Finger Mounting Plate (in-line)
	1F-MP	1 Finger Mounting Plate (in-line)

HINGES/CONNECTORS

 	PRODUCT NUMBER	DESCRIPTION
	GAFFNEY	Gaffney Joint/Hinge (Colt)
	BOA	BOA Fastener for M-Fingers
	DCM 1	1-Finger Direct Connect Mechanism (DCM) Black
	DCM 2	2-Finger Direct Connect Mechanism (DCM) Black
	DCM 3	3-Finger Direct Connect Mechanism (DCM) Black
	DCM 4	4-Finger Direct Connect Mechanism (DCM) Black

GRIP TIPS

	PRODUCT NUMBER	DESCRIPTION
	GRIP TIP	Silicone protective covering , 5 pack

REPLACEMENT PARTS

PRODUCT NUMBER	DESCRIPTION
PMF-BASE-BLK	Partial M-Finger base (Black)
PMF-BASE-BRN	Partial M-Finger base (Brown)
PMF-BASE-TAN	Partial M-Finger base (Tan)
TB-TIP-BLK	Thumb Tip (Black)
TB-TIP-BRN	Thumb Tip (Brown)
TB-TIP-TAN	Thumb Tip (Tan)
TB-BASE-BLK	Base of thumb (Black)
TB-BASE-BRN	Base of thumb (Brown)
TB-BASE-TAN	Base of thumb (Tan)
SM-MJ-BLK	Small Middle joint (Black)
SM-MJ-BRN	Small Middle Joint (Brown)
SM-MJ-TAN	Small Middle Joint (Tan)
MED-MJ-BLK	Medium Middle joint (Black)
MED-MJ-BRN	Medium Middle Joint (Brown)
MED-MJ-TAN	Medium Middle Joint (Tan)
LG-MJ-BLK	Large Middle joint (Black)
LG-MJ-BRN	Large Middle Joint (Brown)
LG-MJ-TAN	Large Middle Joint (Tan)
FB-BLK	Finger Base Joint (Black)
FB-BRN	Finger Base Joint (Brown)
FB-TAN	Finger Base Joint (Tan)
SM-FT-BLK	Small Finger Tip (Black)
SM-FT-BRN	Small Finger Tip (Brown)
SM-FT-TAN	Small Finger Tip (Tan)
MED-FT-BLK	Medium Finger Tip (Black)
MED-FT-BRN	Medium Finger Tip (Brown)
MED-FT-TAN	Medium Finger Tip (Tan)
LG-FT-BLK	Large Finger Tip (Black)
LG-FT-BRN	Large Finger Tip (Brown)
LG-FT-TAN	Large Finger Tip (Tan)



Partial Hand SOLUTIONS, LLC

Distributed by:



college park

800.728.7950 | 586.294.7950

Learn more at college-park.com/partial-hand-solutions

